


Nanodyne Replacement Illuminator for Nikon E600 - Included Items

Additional Items Included But Not Shown:

- PN 10456 Hex Key 1.5mm (for pot knob)
- PN 10471 Hex Key 2mm x 9 inch T-Handle (for securing illuminator to scope)

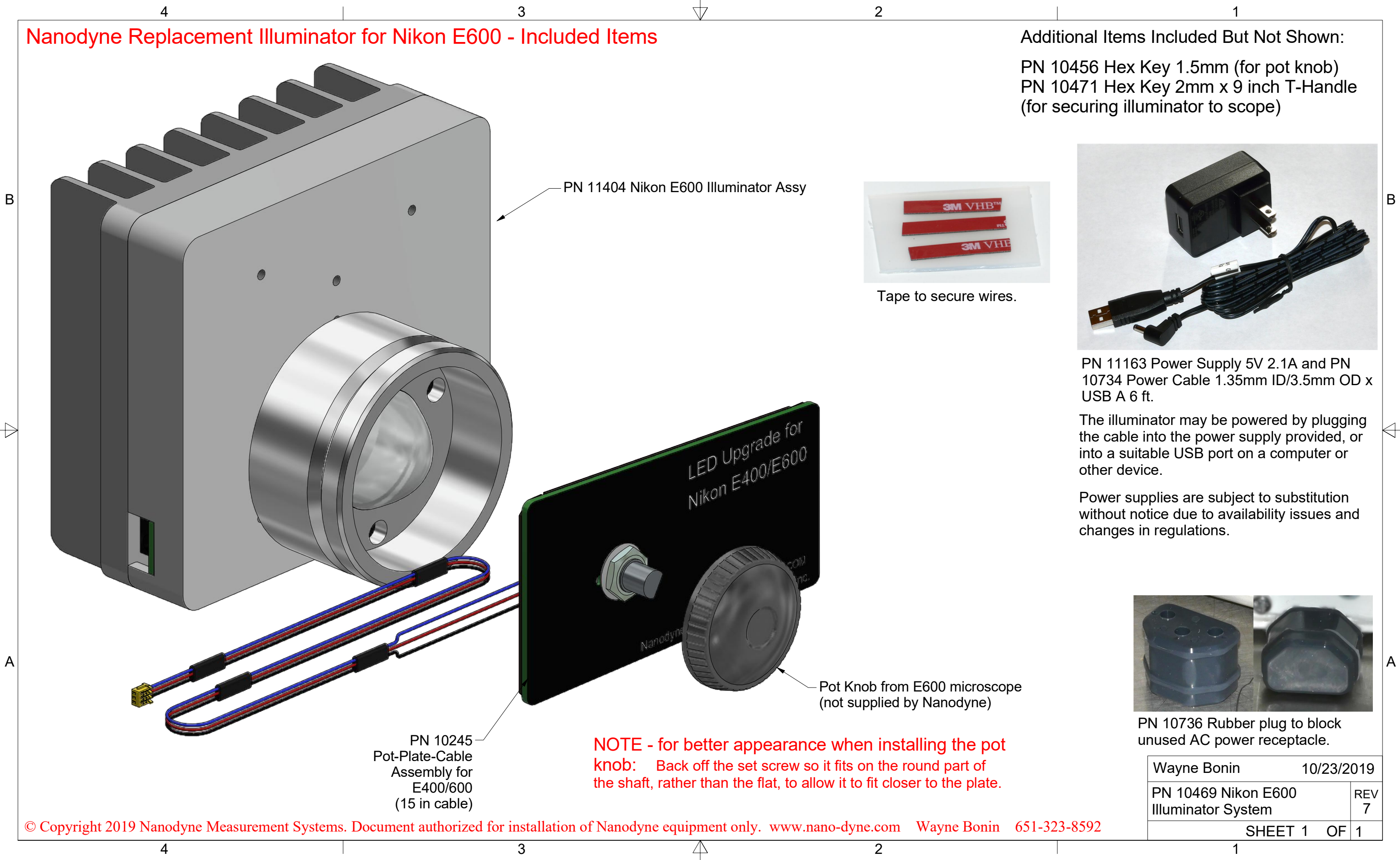

PN 11404 Nikon E600 Illuminator Assy



Tape to secure wires.



PN 11163 Power Supply 5V 2.1A and PN 10734 Power Cable 1.35mm ID/3.5mm OD x USB A 6 ft.

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.

Power supplies are subject to substitution without notice due to availability issues and changes in regulations.

PN 10245 Pot-Plate-Cable Assembly for E400/600 (15 in cable)

NOTE - for better appearance when installing the pot knob: Back off the set screw so it fits on the round part of the shaft, rather than the flat, to allow it to fit closer to the plate.

Pot Knob from E600 microscope (not supplied by Nanodyne)



PN 10736 Rubber plug to block unused AC power receptacle.

Wayne Bonin	10/23/2019
PN 10469 Nikon E600 Illuminator System	REV 7
SHEET 1	OF 1