


# Nanodyne Replacement Illuminator for Nikon 80i - Included Items

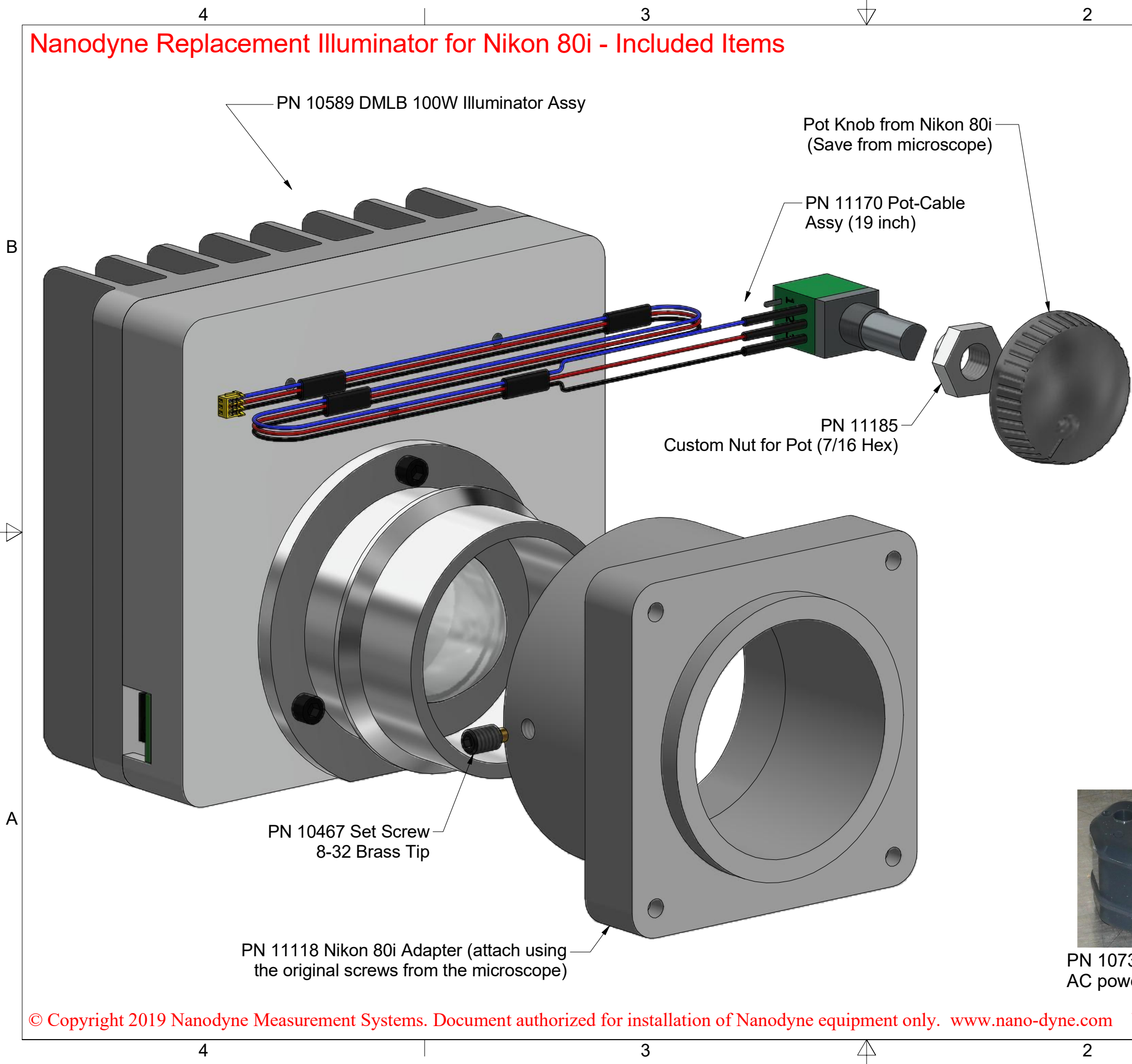

## Additional Items Included But Not Shown:

- PN 10456 Hex Key 1.5mm (for pot knob)
- PN 10457 Hex Key 2.5mm Short Arm L
- PN 10337 Hex Key 2mm Short Arm (for set screw)

## Needed but not included:

- 12mm socket to remove original pot
- 7/16 socket to install new pot

**IMPORTANT NOTE** - If the hole for the potentiometer is less than 9mm diameter you need to drill it out with a 9mm or 23/64 drill to fit the custom nut.

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.



PN 11163 Power Supply 5V 2.1A and PN 10734 Power Cable 1.35mm ID/3.5mm OD x USB A 6 ft.

Power supplies are subject to substitution without notice due to availability issues and changes in regulations.



PN 10736 Rubber plug to block unused AC power receptacle.



3M VHB tape to secure pot cable

Wayne Bonin	12/8/2019
PN 11121 Nikon 80i Illuminator System	REV 3
SHEET 1 OF 1	