


# Nanodyne Replacement Illuminator for Nikon Diaphot 200 & 300 100W - Included Items

PN 11145 T2PIHP3 Illuminator Assy  
w Offset Lens - Side Adjust

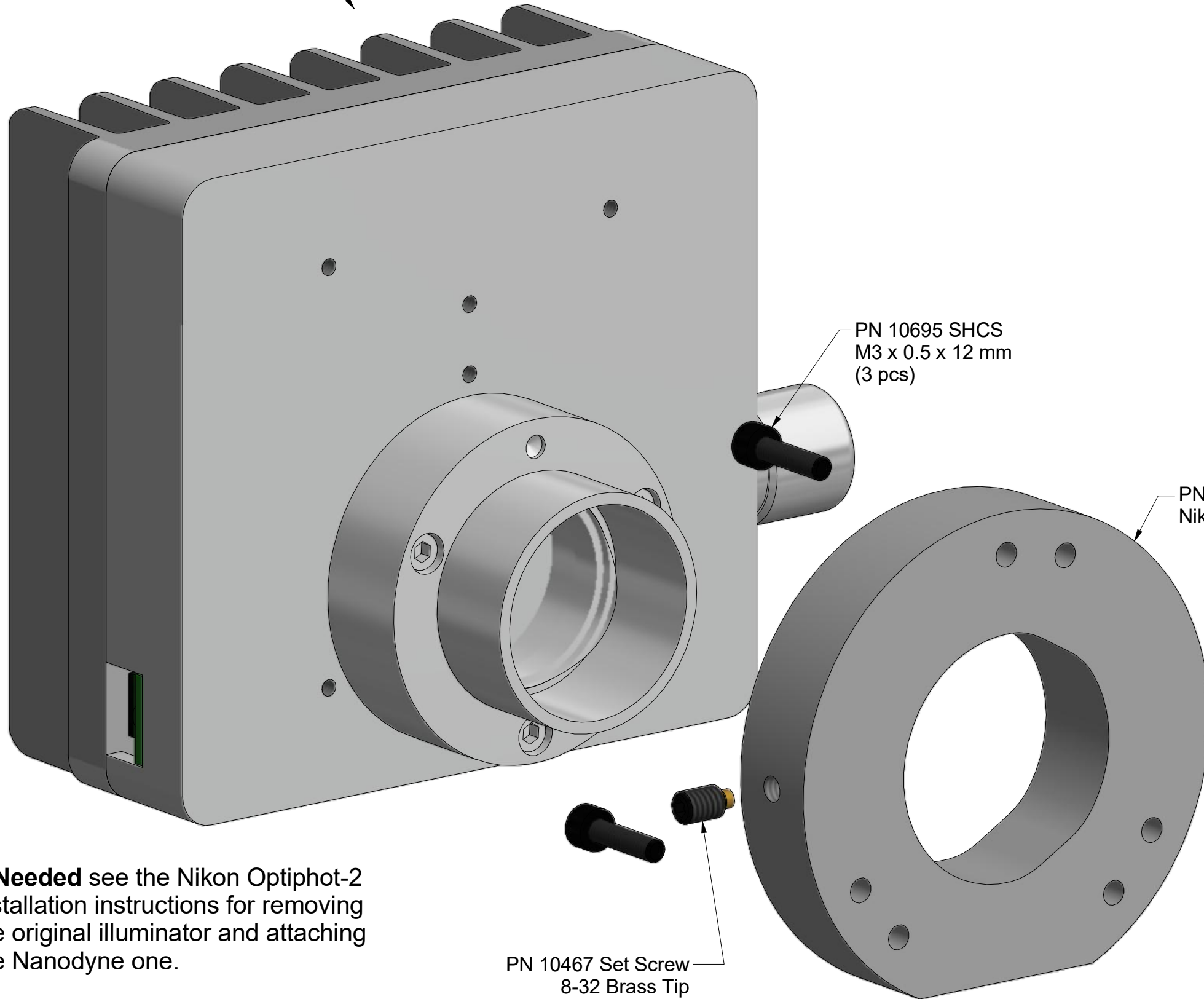


**If Needed** see the Nikon Optiphot-2 installation instructions for removing the original illuminator and attaching the Nanodyne one.

PN 10467 Set Screw  
8-32 Brass Tip

## Additional Items Included But Not Shown:

- PN 10490 Hex Key 5/64 inch (for securing illuminator to adapter)
- PN 10566 Hex Key 2.5mm long arm ball end (for M3 SHCS to attach adapter to microscope)

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.

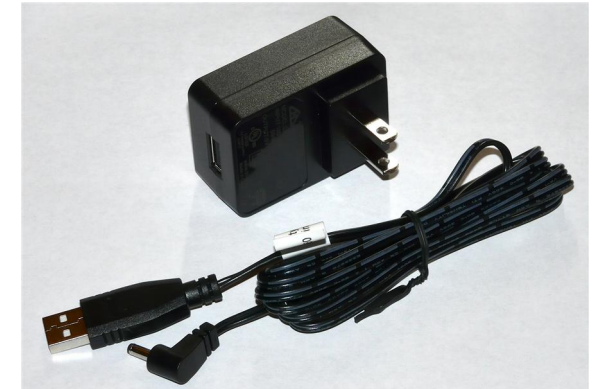

PN 11163 Power Supply 5V 2.1A and  
PN 10734 Power Cable 1.35mm  
ID/3.5mm OD x USB A 6 ft.

Power supplies are subject to substitution without notice due to availability issues and changes in regulations.



PN 10736 Rubber plug to block  
unused AC power receptacle.



3M VHB tape to secure pot cable

Han-Seung Yang	11/5/2020
PN 11226 Nikon Diaphot 200 & 300 - 100W Illuminator System - Side adjust	REV 2 SHEET 1 OF 1