


Nanodyne Replacement Illuminator for Nikon E800 - Included Items

Additional Items Included But Not Shown:
 PN 10456 Hex Key 1.5mm (for pot knob)
 PN 10471 Hex Key 2mm x 9 inch T-Handle
 (for securing illuminator to scope)

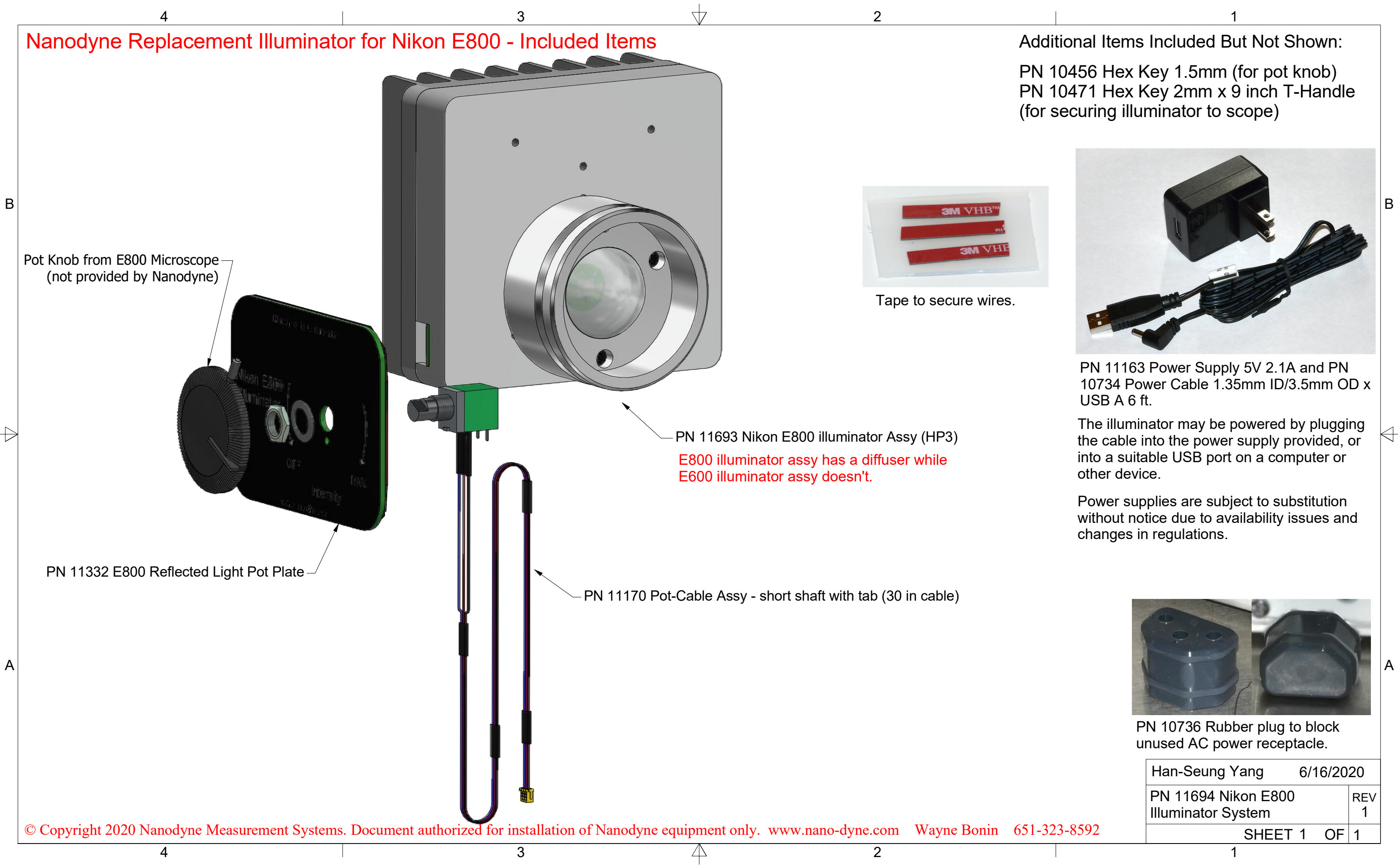

Pot Knob from E800 Microscope
 (not provided by Nanodyne)

PN 11332 E800 Reflected Light Pot Plate

PN 11693 Nikon E800 illuminator Assy (HP3)

**E800 illuminator assy has a diffuser while
 E600 illuminator assy doesn't.**

PN 11170 Pot-Cable Assy - short shaft with tab (30 in cable)



Tape to secure wires.



PN 11163 Power Supply 5V 2.1A and PN 10734 Power Cable 1.35mm ID/3.5mm OD x USB A 6 ft.

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.

Power supplies are subject to substitution without notice due to availability issues and changes in regulations.



PN 10736 Rubber plug to block unused AC power receptacle.

Han-Seung Yang		6/16/2020
PN 11694 Nikon E800 Illuminator System		REV 1
SHEET 1		OF 1