


Nanodyne Replacement Illuminator for Olympus CX31 / 41 Microscope: Included Items

Also Included but not shown:

- PN 10337 Hex Key 2mm Short Arm (for M3 BHCS)
- PN 10457 Hex Key 2.5mm Short Arm (for M3 SHCS)

Tools not included but required:

- Phillips Screwdriver
- 11 mm socket
- Misc tools for removing panels and original illuminator from microscope.

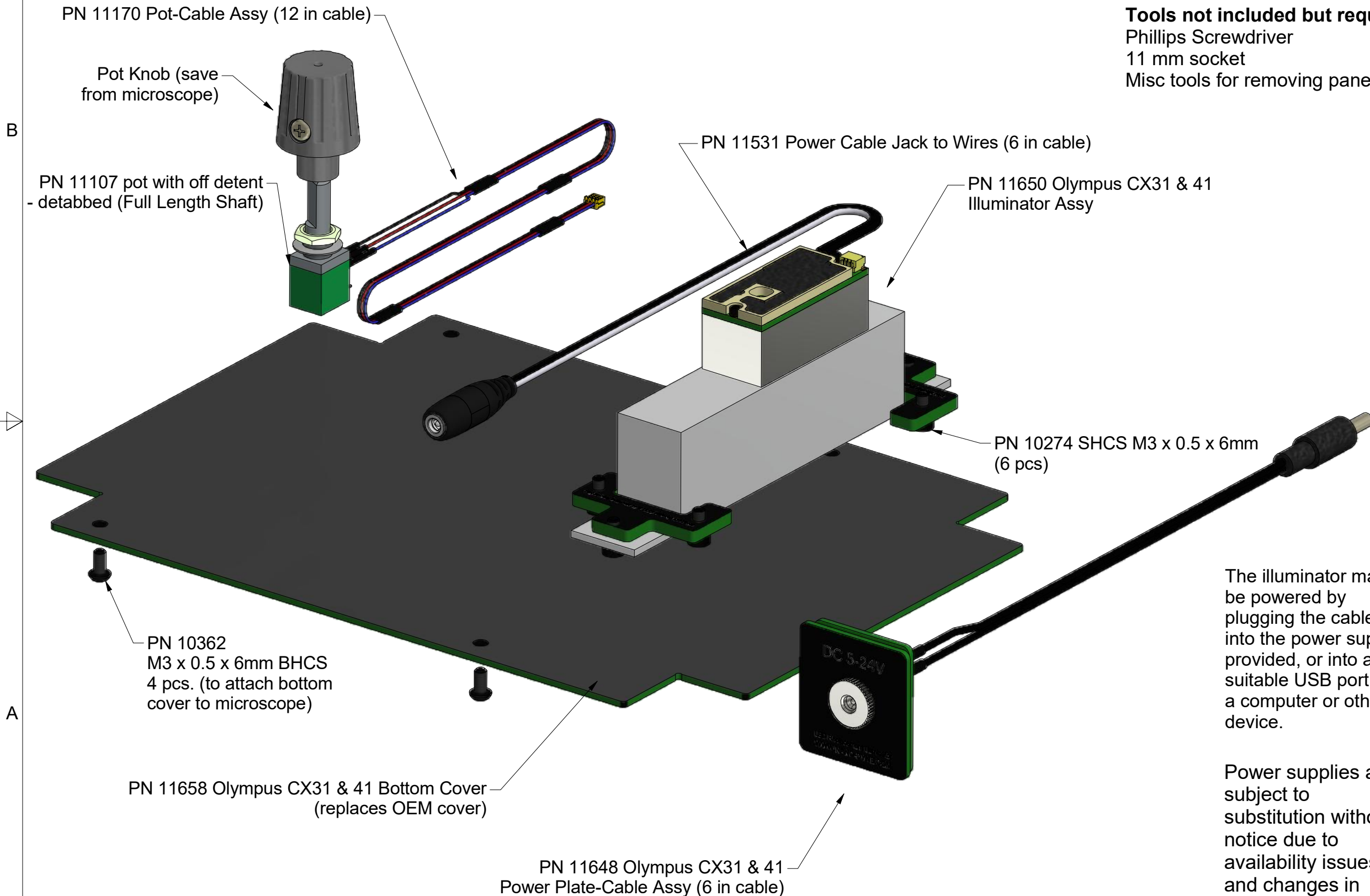

3M VHB Tape to secure wires.

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.



Power supplies are subject to substitution without notice due to availability issues and changes in regulations.

PN 10733 Power Supply & PN 10734 Cable Assy 1.35mm ID x 3.5mm OD RA plug to USB A 6 ft.

Han-Seung Yang	2/4/2020	
PN 11651 Olympus CX31 & 41 Illuminator System		REV 1
SHEET 1 OF 1		