


Nanodyne Replacement Illuminator for Olympus BX51 Bottom: Included Items

Additional Items Included But Not Shown:

PN 10457 Hex Key 2.5mm Short Arm (for M3 screws)
PN 10455 Hex Key 3/32 Ball End L-Key Stubby (for 4-40 screws)
PN 10275 SHCS M3 x 0.5 x 8mm (3pcs) - in case your condenser lens is NOT as pictured on sheet 3.
PN 10697 Cable Ties.

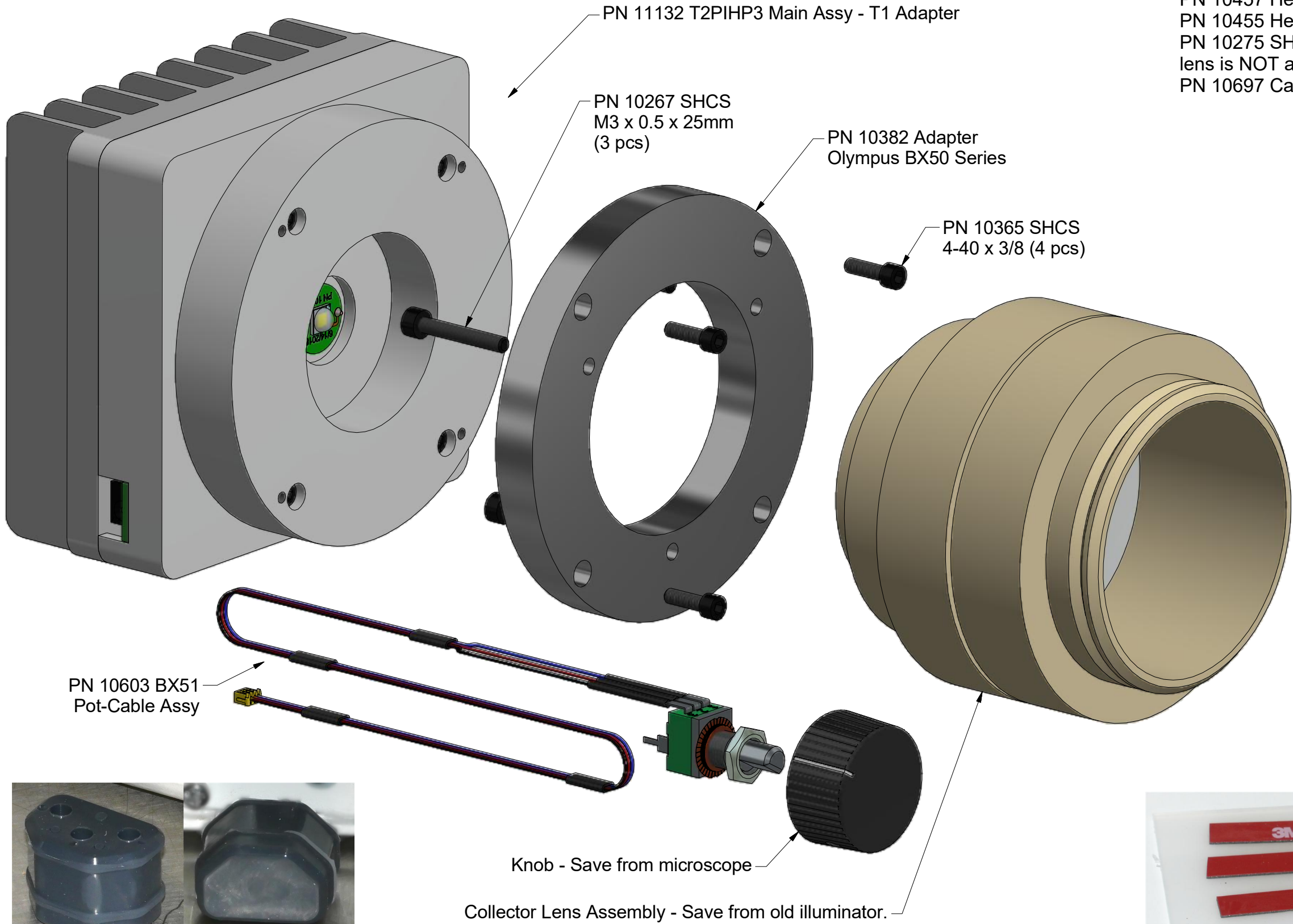

PN 11163 Power Supply 5V 2.1A and PN 10734 Power Cable 1.35mm ID/3.5mm OD x USB A 6 ft.

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.



PN 10736 Rubber plug to block unused AC power recepticle.



3M VHB tape to secure pot cable

Power supplies are subject to substitution without notice due to availability issues and changes in regulations.

Wayne Bonin		2/7/2019	
PN 10413 Olympus BX51 Bottom Illuminator System		REV	4
SHEET 1		OF	1