

4

3

2

1


Included Items Not Shown:

- PN 10457 2.5mm Hex Key (for 3 mm screws to attach OEM collector lens)
- PN 10456 Hex Key 1.5mm (to attach new slide pot to microscope)

The illuminator may be powered by plugging the cable into the power supply provided, or into a suitable USB port on a computer or other device.



PN 11163 Power Supply 5V 2.1A and
PN 10734 Power Cable 1.35mm
ID/3.5mm OD x USB A 6 ft.

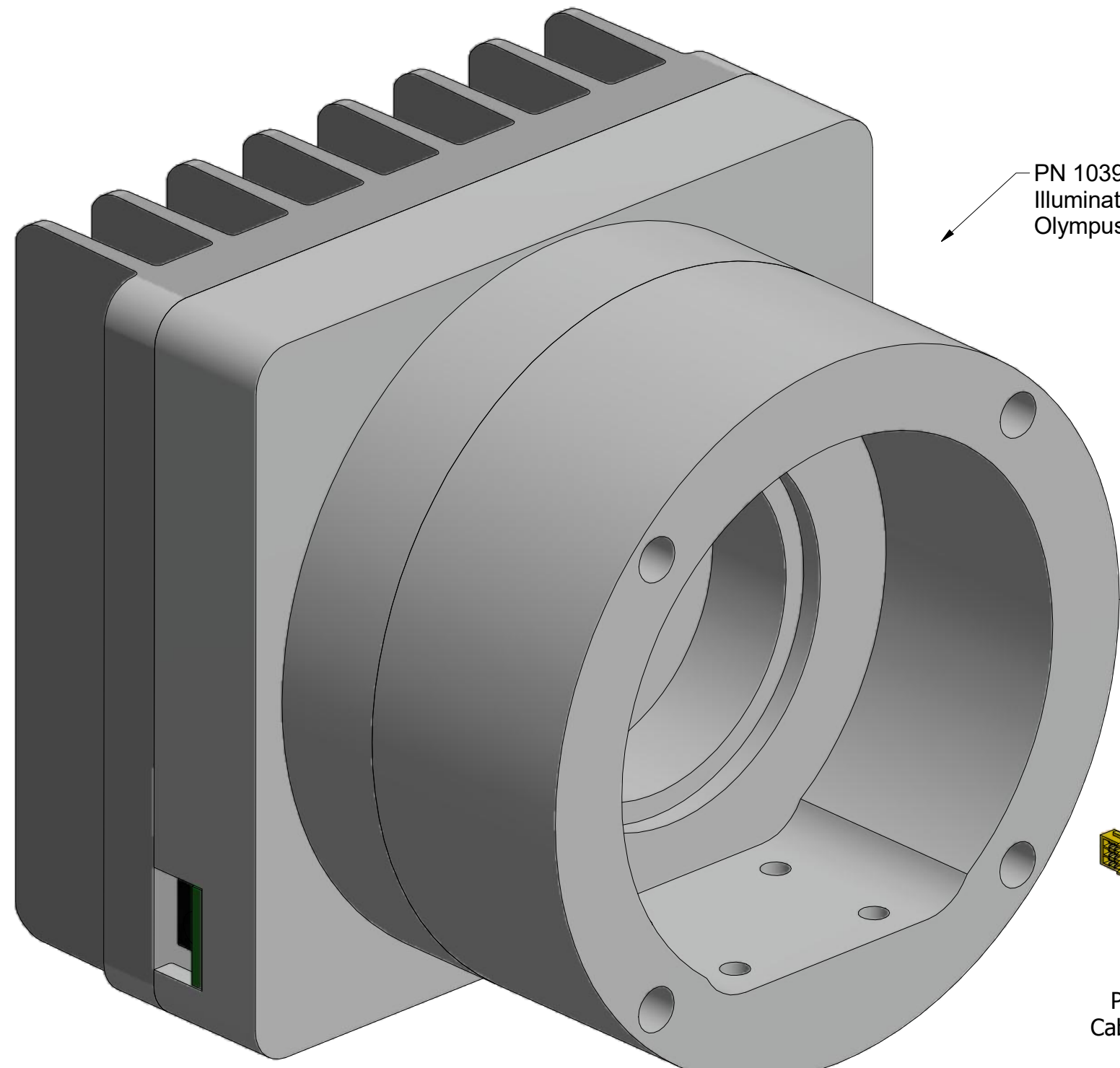

PN 10393 HP3
Illuminator Assy
Olympus BHS

Pot Knob not supplied
(save from microscope)

PN 10276
Slide Pot
Assembly

PN 11002 Pot
Cable Assembly
(21 inch)



Tape to secure wires.



PN 10736 Rubber plug to block unused
AC power receptacle.



PN 10275 SHCS M3 x 0.5 x 8mm (4 pcs)
(to attach OEM collector lens to illuminator)

Power supplies are subject to
substitution without notice due to
availability issues and changes
in regulations.

Wayne Bonin	2/4/2019		
PN 10476 Olympus BHS Illuminator System		REV	5
SHEET 1	OF	1	